


Automated Tissue stainer MAS

Linear Multistainer –
High efficiency for modern labs

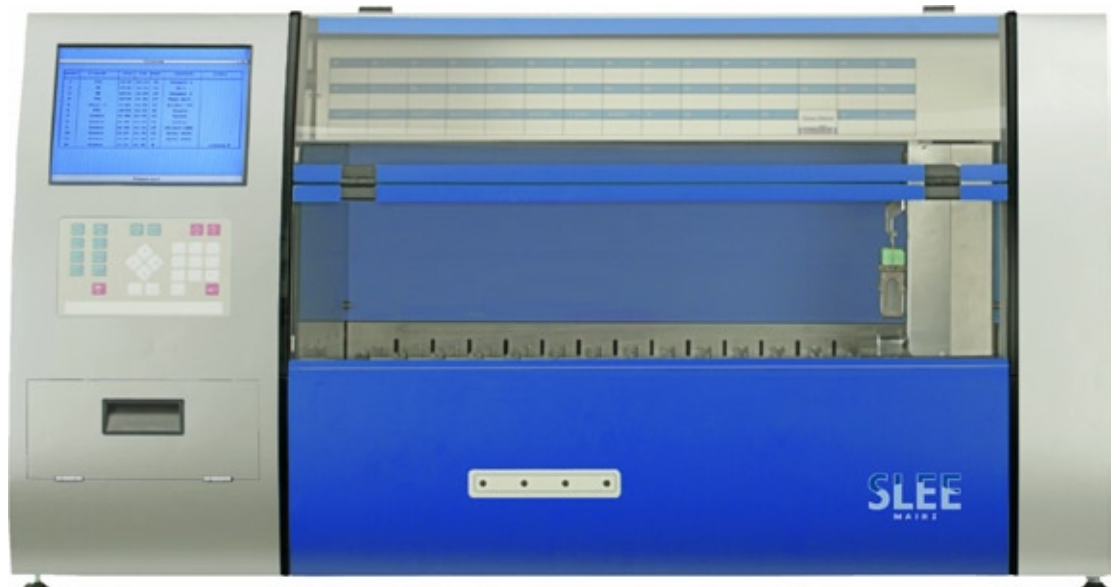

Microtomes

Cryostats

Paraffin

Embedding Centers

Tissue processors

Stainers

MSM

► **MAS**

Cover slipper

Knife sharpener

Floating bath

Slide warmer

Automated Tissue stainer MAS

Linear Multistainer –

High efficiency for modern labs

SLEE
M A I N Z

The tissue stainer MAS combines flexibility with speed and opens new dimensions in the field of automation. It is possible to work with different staining protocols at the same time due to the large number of stations (44). Use of up to 12 baskets (each 30 slides), LCD screen shows all programmes and information, coloured display, safety system for reagents, water safety function, integrated fan for user safety

- ▶ 44 stations
 - Consisting of 4 loading stations
 - 4 unloading stations
 - 6 water baths for rinsing water
 - 30 staining stations with agitation
- ▶ 2 drying stations optional (in combination with 4 waterbaths)
- ▶ Active carbon filter or the tubing for a sucking system is best for user safety.
- ▶ Universal power input 100 – 240 V, 50/60 Hz
- ▶ Optional:
 - 2 drying stations (in combination with 4 water baths)



Stainer MAS
in combination with
Cover Slipper MCS



Stainer MSM for universal use

Technical Data

| | |
|----------------------------------|--------|
| Max. No. of staining protocols | 20 |
| Steps for each protocol | 40 |
| Volume per bath | 360 ml |
| With option for volume reduction | 260 ml |

| Weight and Dimensions | Depth x Width x Height |
|---------------------------|------------------------|
| Dimensions (in mm) | 670 x 1210 x 640 |
| Dimensions packed (in mm) | 880 x 1330 x 770 |

| | |
|---------------|--------|
| Weight net | 115 kg |
| Weight packed | 146 kg |

| | | |
|-----------------------|-------------------|---------------|
| Electrical parameters | Art.-No. 14000500 | 230 V / 50 Hz |
| | Art.-No. 14000501 | 230 V / 60 Hz |
| | Art.-No. 14000502 | 115 V / 60 Hz |
| | Art.-No. 14000503 | 100 V / 50 Hz |
| | Art.-No. 14000504 | 100 V / 60 Hz |

Please ask for other units for Histopathology and consumables as well as laminar flow equipment.

(Specifications are subject to change without notice)

SLEE MEDICAL GMBH

P. O. Box 10 05 07
D-55129 Mainz · Germany

Fon +49 (0) 6131/95871-0
Fax +49 (0) 6131/95871-0

mail@sleemainz.com
www.sleemainz.com